


**Product Name :**  
Gillmore Apparatus

**Product Code :**  
TEST5054RFQLAB-0004

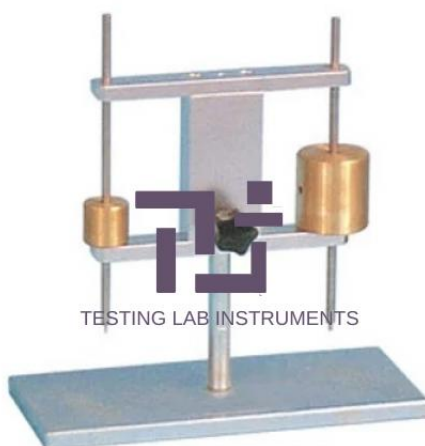

**Description :**

Gillmore Apparatus

**Technical Specification :**

For determination of setting time of hydraulic cement.

The apparatus comprises the following :-

One base with a central vertical rod having two horizontal cross-arms. The lower cross arm is adjustable in height.

Initial Needle 2.12 .mm dia, 113.4 g weight.

Final Needle 1.06 mm dia, 453.6 g weight.

One glass plate.

We are a leading manufacturers of Gillmore Apparatus testing equipment for a range of industries, including materials science testing lab. We offer full range of Gillmore Apparatus testing equipment, including Gillmore Apparatus testing kit, Gillmore Apparatus testing machine, Gillmore Apparatus testing tools, Gillmore Apparatus test equipment, Gillmore Apparatus test machines, Gillmore Apparatus for Lab Testing Instruments. We are suppliers of materials testing tools and equipment for the construction industry. Our mission is to provide reliable, accurate, and easy-to-use Gillmore Apparatus testing equipment that helps customers build better, safer, and more sustainable infrastructure.

---

## Testing lab instruments India

