


**Product Name :**  
Stiff Frame

**Product Code :**  
TEST5045RFQLAB-0001

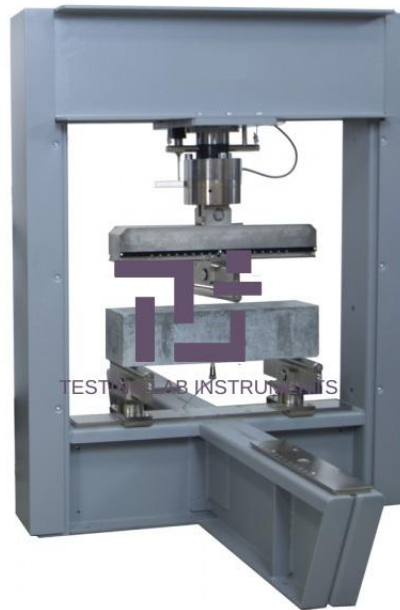

**Description :**

Stiff Frame for Testing Lab Equipment Tender Suppliers India

**Technical Specification :**

The system is a totally welded structure, which gives a tremendous amount of stability to the system. It has an equal area hydraulic cylinder attached to the base and an oil filled spherical seating is screwed to the top with Solid platen having hardness greater than The horizontal and vertical clearances are enough to test standard rock specimen under triaxial test. It can accommodate a triaxial cell with a sample size up to 90 mm size. Front cover is provided, made of acrylic sheet, as a protection to the operator while at the same time giving an unobstructive view of the specimen under test.

Loading Capacity of the unit  
(Static & Dynamic) : 1000/2000kN  
Total travel of the ram : 50mm  
Clearance – Horizontal : 450mm  
Vertical : 820mm  
Platen Size : 250mm

We are a leading manufacturers of Stiff Frame testing equipment for a range of industries, including materials science testing lab. We offer full range of Stiff Frame testing equipment, including Stiff Frame testing kit, Stiff Frame testing machine, Stiff Frame testing tools, Stiff Frame test equipment, Stiff Frame test machines, Stiff Frame for Lab Testing Instruments. We are suppliers of materials testing tools and equipment for the construction industry. Our mission is to provide reliable, accurate, and easy-to-use Stiff Frame testing equipment that helps customers build better, safer, and more sustainable infrastructure.



# Testing lab instruments India